



Sovereign SM6 High Care Ventraflow Pressure Care Mental Health Mattress

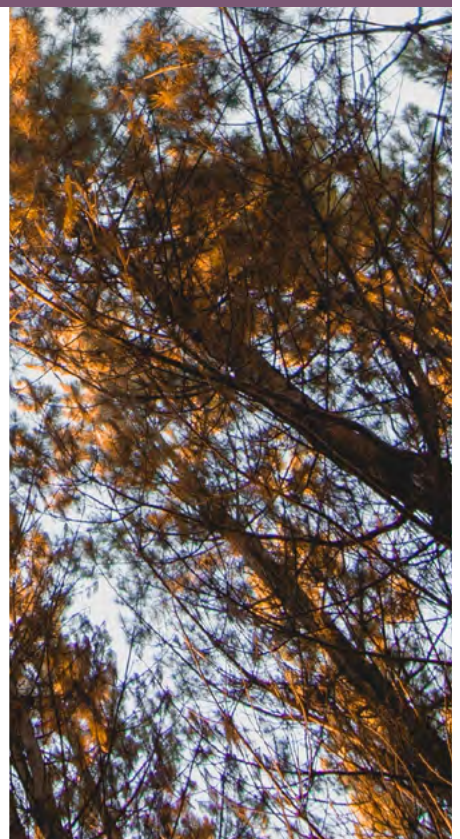


MEDICAL MATTRESS EXPERTS

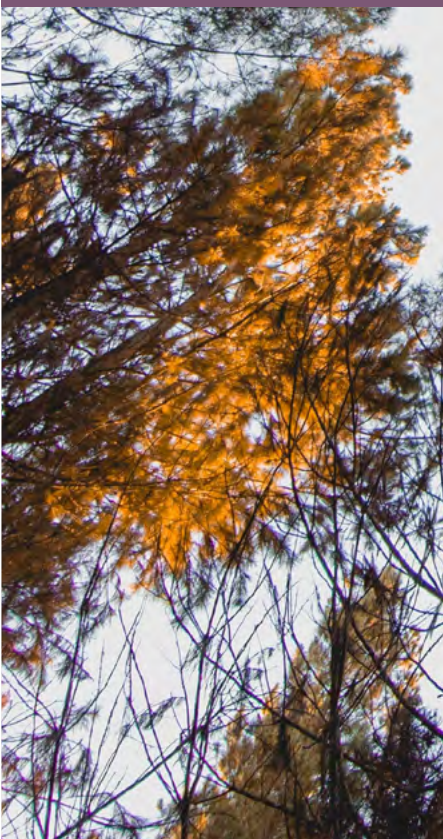

The Sovereign range is the Ultimate blend of longevity, patient care and comfort. Versatility, safety and value is delivered Through its unique and extensive features. Gentle pressure dispersion and redistribution through unique foam technology is integrated together with Foundational and Structural Support for Safer Patient Transfer and Repositioning. This combines to deliver the most advanced, high-performance static Support surface available on today's medical mattress market.


The Sovereign Mental Health range is the benchmark in high risk mental health mattress by combining this high level of pressure care with mental health standards. All covers are fully sealed and enclosed covers to protect the foam core, and all foams used meet fire retardant standards for compliance within Mental Healthcare.

Features

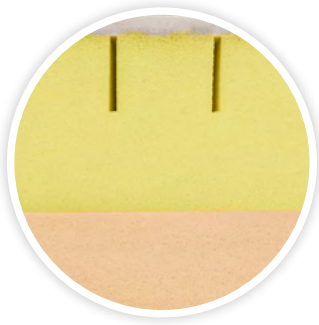

HIGH PERFORMANCE, MENTAL HEALTH SPECIFICATION CORE

- ✓ High specification foam for longevity and performance
- ✓ Provides graded postural support and patient comfort
- ✓ Anti-microbial treatment applied to foam maximises longevity and core micro-climate
- ✓ Conforms to EPUAP, NPUAP & PPPIA Guidelines for quality and comfort
- ✓ Fire retardant - meets BS5852 source 5 specification
- ✓ Profiles & conforms well to Acute Care and Mental Health Beds

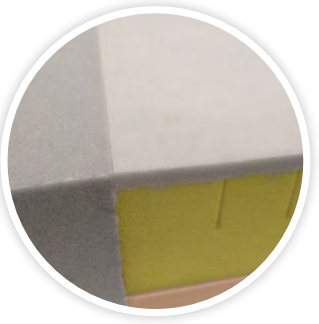


AIROWORM RESPONSE IMMERSION INTERFACE LAYER

- ✓ A particularly responsive Ultra-High Performance elastic foam designed to optimise patient mobility, function and enhance repositioning
- ✓ Open-cell structure provides low-shear surface, good pressure redistribution, effective immersion, enhanced breathability and superior heat dissipation
- ✓ Alleviates localised pressure points
- ✓ Enhances contact with patient for distribution of interface pressure
- ✓ High temperature stability to deliver consistent immersion, performance and comfort
- ✓ Provides exceptional comfort
- ✓ Greatly reduces risk of impression marks on patient skin
- ✓ Firmness: Medium



5-ZONED SECONDARY IMMERSION REGION

- ✓ Delivers enhanced immersion and support for critical body regions i.e. heel/scapular/sacrum
- ✓ Graded support and patient comfort in supine, prone, lateral and fowler positions

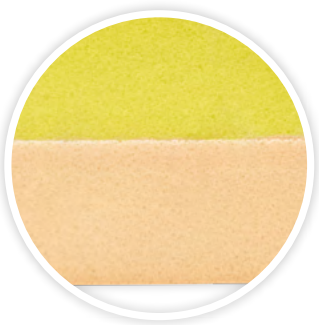

PREMIUM POSTURAL SUPPORT MID-CORE REGION

- ✓ Provides exceptional comfort and postural support for repositioning and patient care
- ✓ Provides exceptional conformity and profiling to medical bed surfaces



PREMIUM FIRM FOUNDATION LAYER

- ✓ Provides protection against 'bottoming-out' during transfer, repositioning and whilst in the inclined/fowler position
- ✓ Provides exceptional comfort, support for repositioning and patient care




HINGED STRENGTHENED SIDES

- ✓ Increase lateral stability of support surface
- ✓ Improves patient safety and encourages central patient positioning
- ✓ Make patient transfers easier
- ✓ Unique hinge profile reduces foam fatigue and allow mattress to easily conform to varied bed positions
- ✓ Enables ease of primary patient care

Sovereign SM6

Sovereign SM6 High Care Ventraflow Pressure Care Mental Health Mattress

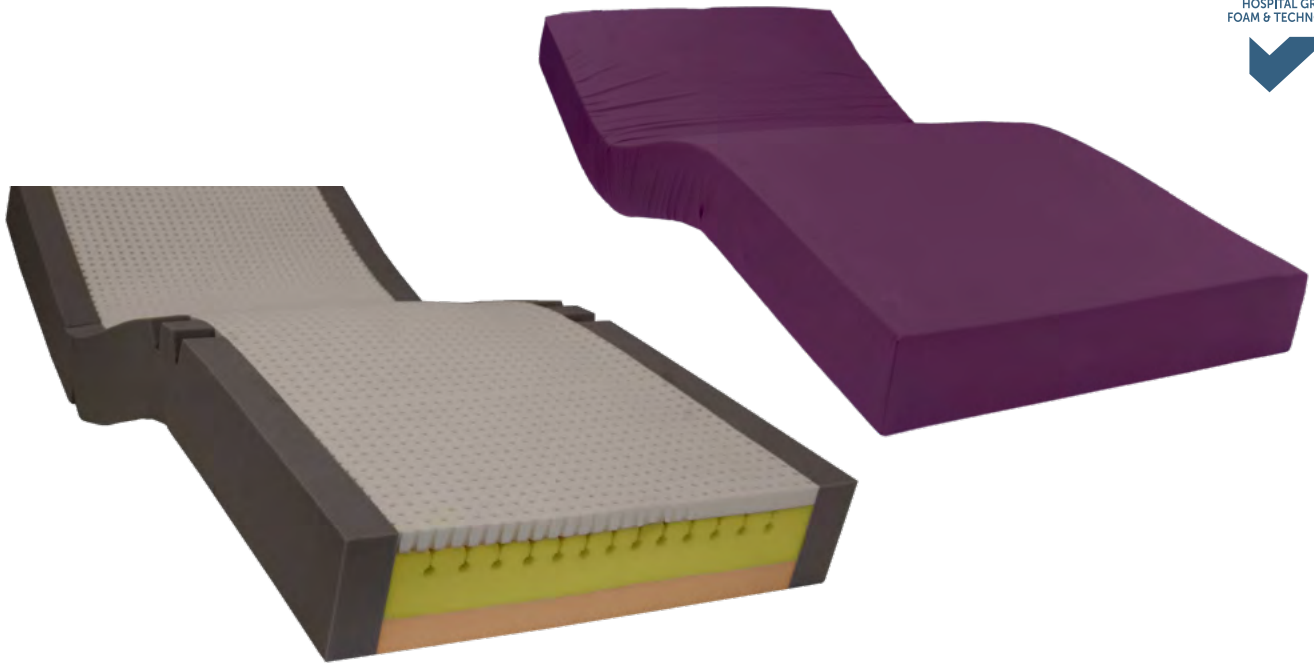

- SHEAR & FRICTION:** ✓ ✓ ✓ ✓ ✓
- REPOSITIONING:** ✓ ✓ ✓ ✓ ✓
- INFECTION CONTROL:** ✓ ✓ ✓ ✓ ✓
- FIRMNESS:** ✓ ✓ ✓ ✓ ✓

- MOISTURE & MICRO-CLIMATE:** ✓ ✓ ✓ ✓ ✓
- PERFORMANCE & LONGEVITY:** ✓ ✓ ✓ ✓ ✓
- MOBILITY & FUNCTION:** ✓ ✓ ✓ ✓ ✓
- PRESSURE CARE RATING:** **MEDIUM-HIGH RISK**

/ Sovereign SM6

CODE	SIZE (MM)	
SM6S1000	1980x880x150	Single
SM6SL1000	2030x880x150	
SM6KS1000	1980x1050x150	King Single
SM6KSL1000	2030x1050x150	
SM6EW1000	1980x1150x150	Extra-Wide
SM6EWL1000	2030x1150x150	
SM6D1000	1980x1350x150	Double
SM6DL1000	2030x1350x150	
SM6Q1000	1980x1530x150	Queen
SM6QL1000	2030x1530x150	
SM6K1000	1980x1850x150	King
SM6KL1000	2030x1850x150	
	Customisable	CAN BE MADE TO SUIT ANY BED
THERAPEUTIC WEIGHT RANGE*	165KG	
APPLICATIONS	<input checked="" type="checkbox"/> MENTAL HEALTH WARDS <input checked="" type="checkbox"/> HOSPITAL WARDS <input checked="" type="checkbox"/> HOME CARE	
WARRANTY	5-YEAR FOAM CORE WARRANTY 4-YEAR COVER WARRANTY (12-MONTH WARRANTY IN MENTAL HEALTH APPLICATION)	
STANDARD COVER	FULLY SEALED & ENCLOSED PREMIFLEX CLASSIC	