



forté
HEALTHCARE


AIR MATTRESS SURROUNDS

AMS1000 Forté Air System Surrounds

Excellent Value, Enhances Functionality

EVERY PRODUCT
CUSTOM MADE

Forté Air Mattress Surrounds provide additional width and/or length to pressure relieving air support surfaces enabling the Air Support Surface to adapt to almost any bed base or configuration. High specification foam enhances the stability of the support surface and provide additional patient and caregiver assistance for positioning and transfers. Covering materials and construction methods match the clinical performance characteristics of the Air Support Surface to ensure that pressure relief and patient skin protection is not compromised. Speak with Forté Healthcare or your distributor today to enquire further.

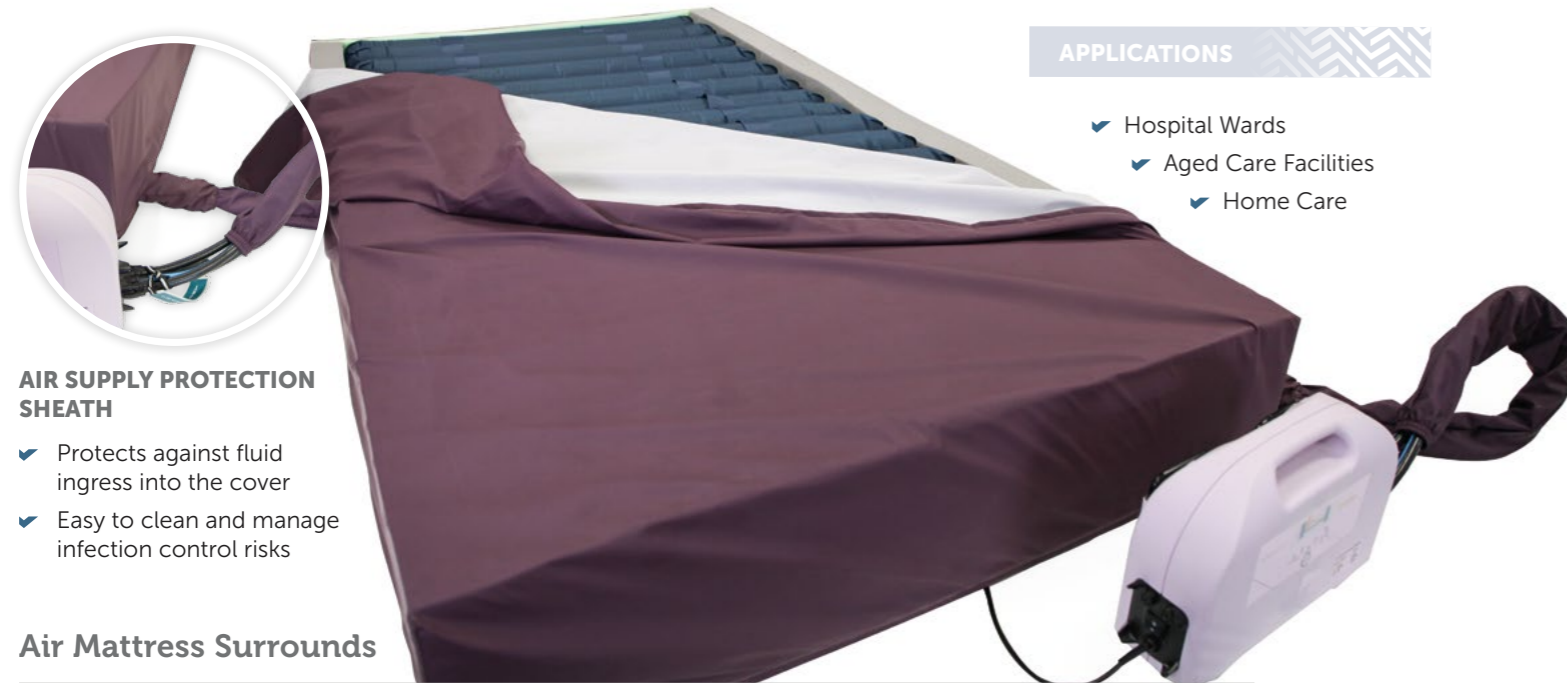

APPLICATIONS

- ✓ Hospital Wards
- ✓ Aged Care Facilities
- ✓ Home Care

AIR SUPPLY PROTECTION SHEATH

- ✓ Protects against fluid ingress into the cover
- ✓ Easy to clean and manage infection control risks

Air Mattress Surrounds

BRAND	CODE	OVERALL SIZE	
Custom	AMS1002		
Harvest Cavalier	AMS1003		
Arjo Nimbus	AMS1004		
Novus Premium 8	AMS1005		
Novus Premium 10	AMS1006	King Single	2000x1050
Talley Quattro Plus	AMS1007		
Talley Quattro Acute	AMS1008		
Talley Pulsair Choice	AMS1009		
Easy Air 8000	AMS1010		
	AMS1011		
Custom	AMS1012		
Harvest Cavalier	AMS1013		
Arjo Nimbus	AMS1014		
Novus Premium 8	AMS1015		
Novus Premium 10	AMS1016	Extra-Wide	2000x1150x200
Talley Quattro Plus	AMS1017		
Talley Quattro Acute	AMS1018		
Talley Pulsair Choice	AMS1019		
Easy Air 8000	AMS1020		
	AMS1021		



SUPPORTIVE FOAM BASE

- ✓ Provides additional protection to air cells
- ✓ Should failure occur within air mattress, foam provides additional protection between patients body and bed base
- ✓ Eliminates the requirement for an additional mattress beneath the surround making the support surface easier to clean



PROFILE STABILISING SIDES

- ✓ Enables ease of primary patient care
- ✓ Facilitates patient transfers
- ✓ Unique hinge points allow mattress to easily conform to varied bed positions
- ✓ Encourages central patient positioning

4-WAY STRETCH, VAPOUR PERMEABLE, PRESSURE CARE UPPER

- ✓ Breathable
- ✓ Replaces air mattress cover
- ✓ High specification to provide maximum shear relief
- ✓ Available in Premiflex Grade
- ✓ Wipe clean
- ✓ Fire retardant, anti-microbial, resistant to body fluids

HEAVY DUTY PVC BASE

- ✓ Provides enhanced durability
- ✓ Delivers great durability & performance
- ✓ Fire retardant, anti-microbial, resistant to body fluids



360° ZIPPER AND WATERFALL INGRESS PROTECTION FLAP

- ✓ Full removable upper cover enables easy installation of air mattress on-site
- ✓ Easy access for internal mattress adjustment/repair
- ✓ Double coated Premiflex Waterfall flap provides protection against fluid ingress. Backing is sealed for ease of cleaning



COMPATIBLE WITH CPR FUNCTION

- ✓ Custom made to provide access to the CPR function. Note - the access style differs by air mattress brand
- ✓ Cover print ensures visibility for training and in-service purposes

MEDICAL MATTRESS EXPERTS



P. 1300 73 11 63
W. fortehealthcare.com.au
E. sales@fortehealthcare.com.au