




 **forté** Soft Touch Ultra Cover

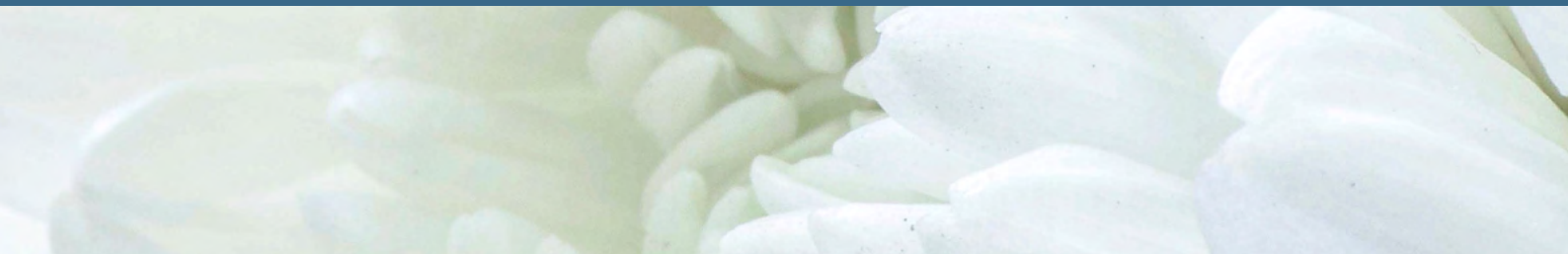

# Soft Touch Ultra Cover

The Soft Touch Ultra cover is a 100% waterproof, durable PVC material manufactured using a casting process. This provides a soft, 'polyurethane type', durable finish. This high performing cover is designed for high demanding patient care scenarios from inpatient mental health to group home care. This cover is designed as a damage resistant, waterproof cover and does offers only limited pressure care.





### SOFT TOUCH ULTRA PVC1150 MATERIAL

- ✓ 100% Waterproof
- ✓ Heavy Duty, Robust & Hardwearing
- ✓ Non abrasive, smooth, no texture
- ✓ Protects foam
- ✓ Superior resistance to cleaning agents
- ✓ Standard Colour: Beige

### CLEANING AND CARE

- ✓ 10,000ppm disinfection chemical suitability
- ✓ Wipe clean
- ✓ Stain Resistant
- ✓ Anti-microbial protected



### COMPLIANCE

- ✓ Fire retardant, meets BS 6807 ignition source 5 and BS 7175 ignition source 5 standards
- ✓ Anti-microbial protected
- ✓ All seams high frequency welded, no exposed sewn seams
- ✓ Meets and Exceeds Facility Requirement
- ✓ Material manufactured in Belgium, to Forté specification

### CONSTRUCTION

- ✓ All seams high frequency sealed & welded
- ✓ 180° / two-sided integrated zipper
- ✓ Waterfall Fluid Penetration Flap
- ✓ TGA Compliant Product Labelling



### INTEGRATED ZIPPER

- ✓ 180° / Two-Sided Zipper
- ✓ Allows easy cover removal
- ✓ Robust and durable
- ✓ High frequency welded attachment for optimal protection against fluid ingress

### WATERFALL ZIPPER PROTECTION FLAP

- ✓ Provides additional protection against fluid ingress through zipper
- ✓ Meets infection control guidelines



### HIGH FREQUENCY SEALED AND WELDED SEAMS

- ✓ Electromagnetically fused seams protecting against risk of penetration of fluid
- ✓ Meets and exceeds infection control requirements
- ✓ Reinforced with additional heavy duty stitching at stress points
- ✓ 100% Sealed to prevent ingress of fluid or infection

# Options



## FULLY SEALED & ENCLOSED

- ✓ No-Zipper
- ✓ All seams welded and 100% waterproof to prevent ingress of fluid and infection (high frequency sealed and welded seams)
- ✓ Ultimate core durability and protection,
- ✓ Used for mental health - protects foam core against willful damage
- ✓ Used in high risk infection control applications
- ✓ Usually manufactured with an integrated pressure relief system



## CARRY HANDLES

- ✓ Fits to side of Mattress Cover for easy access whilst Mattress is on the bed
- ✓ Internal Reinforcement for Longevity
- ✓ Close Mount Design Significantly Reduces Risk of Getting Caught on Other Patient Care Equipment
- ✓ Fully Sealed Non-Absorbent and Wipe-Clean for Infection Control
- ✓ Very Strong and Flexible for Performance and Durability



## CUSTOM PRINTING

- ✓ Print ward location, facility name, asset number or any other details to track and control where your mattress gets used
- ✓ Audit register
- ✓ Print additional care information specific to your facility
- ✓ Custom branding available
- ✓ Full colour, indelible, digital system providing excellent flexibility and versatility



## EVACUATION SYSTEMS

- ✓ Fully Integrated / Built-in Forté Premium Assure Evacuation System
- ✓ Standard Forté evacuation system accessory
- ✓ Contact Forté directly for further details



## CUSTOM ZIPPER

- ✓ up to 360° / Four-Sided Zipper
- ✓ Allows easy cover removal
- ✓ Robust and durable
- ✓ High frequency welded attachment for optimal protection against fluid ingress
- ✓ Other custom zippers available - contact Forté and provide them with the details of what you require. i.e. water resistant, invisible.

# FORTÉ PVC1150 SOFT-TOUCH WATERPROOF FABRIC

## Typical Applications


Medical mattress covers  
 Mental health applications  
 Support cushions  
 Childcare mattresses  
 Hostels and dormitories

## Features & Benefits

Excellent HF Welding Compatibility	Allows covers to be sealed against fluid ingress at seams
100% Waterproof	Allows care staff to 'Wipe-Clean', effective barrier against fluid ingress
Soft Touch Feel	With the durability of a conventional PVC, Soft Touch Ultra has a softer polyurethane type feel for enhanced comfort
Fire Retardant	Meets hospital and medical facility requirements
Excellent Flexibility and Cracking Resistance	Greater durability and increased infection control



Fire retardant

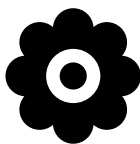

Hypoallergenic



Anti-bacterial



Waterproof



Stain resistant



Soft handle



Anti-Static

Property	Result	Test Method
Width	110, 220cm	ISO 3932
Tearing resistance	Warp: 220N Weft: 200N	AS2001.2.10
Tensile Resistance	Warp: 900N Weft: 1000N	AS2001.2.3
Adhesion	35N/5cm min	AS1441.5
Flame Retardant	Pass	AS1441.13
Flex Cracking	>400,00	AS1441.6
Anti-Microbial	No fungus growth	AS 1157
Washing	75°C	N/A
Dry cleaning	This is Dry cleanable and HF weldable	N/A
<b>Care and Maintenance</b>		
Soft-Touch requires regular care and maintenance. If the following instructions are followed, the end user can expect to maximise the fabrics potential life and maintain a hygienic environment.		
General Cleaning	Regularly wipe with a mild solution of phenolic based cleaner and water. Rinse thoroughly with fresh water and wipe dry	
Spot Cleaning	For difficult spills and stains, prepare a strong solution of phenolic based commercial cleaner. Using a soft bristle brush and solution, clean the affected area. Rinse thoroughly with fresh water and wipe dry.	
Disinfection	It is common to disinfect surfaces in the healthcare environment. Using a water/bleach solution of 10:1, clean the entire mattress. The bleach must not have sodium hypochlorite level in excess of 5%. Disinfectants applied at full concentration will decrease the life of Soft-Touch Materials. Rinse thoroughly with fresh water and wipe dry.	
Notes	<ul style="list-style-type: none"> <li>• Never exceed the cleaning chemical manufacturers concentration recommendations</li> <li>• Do not use abrasive cloths or cleaners</li> <li>• Do not use aggressive cloths or chemicals</li> <li>• Non autoclavable</li> <li>• Contact Forte Healthcare with any further questions: 1300 73 11 63</li> </ul>	