



forté
HEALTHCARE

FINALÉ MEDICAL BED PROTECTORS & ACCESSORIES

BED RAIL PROTECTORS - BOLSTERS - FALL MATS

Finalé BRP1 Wrap Around Bed Rail Protectors



Customisable Wrap Around Bed Rail protectors provide protection to the patient from bed rails. Utilising facility grade fabrics and high-density foams, Forte Bed Rail protectors conform to infection control and cleaning guidelines. Available in a range of fastening methods for simple fitment by care staff. Cover options include Heavy Duty PVC, Pressure Care Material and Mental Health Grade Material.



ADVANTIFLEX PRESSURE CARE COVER

ADVANTIFLEX PU1100 POLYURETHANE PRESSURE CARE MATERIAL:

- ✓ Vapour permeable and breathable
- ✓ Anti-shearing
- ✓ Fire retardant, resistant to body fluids and bacteria
- ✓ Wipe-clean
- ✓ Heavy Duty PVC Inner Material – for maximum durability
- ✓ Optional: heavy duty PVC material for maximum durability

CONSTRUCTION:

- ✓ Closure style: velcro or press - for quick and simple fitment.
- ✓ Robust seams

25MM PADDED PROTECTIVE FOAM CORE

- ✓ High quality, open cell foam
- ✓ Anti-microbial treatment
- ✓ Padded on both sides of rail for maximum patient protection

BRP1 Wrap Around Bed Rail Protectors (pair)

COVER	CLOSURE STYLE	SIZE (MM)	CODE
Advantiflex	Velcro Closure	1500x400	BRP1000
		1800x400	BRP1001
		2000x400	BRP1002
	Press-Stud Closure	1500x400	BRP2000
		1800x400	BRP2001
		2000x400	BRP2002
Heavy Duty PVC	Velcro Closure	1500x400	BRP5000
		1800x400	BRP5001
		2000x400	BRP5002
	Press-Stud Closure	1500x400	BRP5003
		1800x400	BRP5004
		2000x400	BRP5005

APPLICATIONS

- ✓ Hospitals
- ✓ Aged Care
- ✓ Group Homes
- ✓ Home Care

Finalé BRP2 Padded Mesh Bed Rail Protectors



Finalé Padded Mesh Bed Rail Protectors are designed for fitment to full-length bed rails to provide additional patient protection. The integration of a mesh panel provides patient and caregivers with increased visibility and reduces any sense of claustrophobia. Removable and washable, Finalé Padded Mesh Rail protectors are popular with pediatric clients. Custom options include a triple-bar protector and different fabric colours to fit within the care setting.



ADVANTIFLEX PRESSURE CARE COVER

ADVANTIFLEX PU1100 POLYURETHANE PRESSURE CARE MATERIAL:

- ✓ Vapour permeable and breathable
- ✓ Anti-shearing
- ✓ Fire retardant, resistant to body fluids and bacteria
- ✓ Wipe-clean
- ✓ Optional: heavy duty PVC material for maximum durability

MESH INSERT PANEL:

- ✓ Enhances patient and caregiver visibility
- ✓ No compromise of rail functionality and folding
- ✓ Reduces any sense of claustrophobia

CONSTRUCTION:

- ✓ Closure style: Heavy Duty Zipper fitment - Robust & Secure
- ✓ Robust seams

20MM PADDED PROTECTIVE CORE

- ✓ High quality, open cell foam
- ✓ Anti-microbial treatment
- ✓ Padded on both sides of rail for maximum patient protection

BRP2 Padded Mesh Wrap Around Bed Rail Protectors

COVER	SIZE (MM)	CODE
Advantiflex	1500x400	BRP21000
	1800x400	BRP21001
	2000x400	BRP21002
Heavy Duty PVC	1500x400	BRP25000
	1800x400	BRP25001
	2000x400	BRP25002

APPLICATIONS

- ✓ Hospitals
- ✓ Aged Care
- ✓ Group Homes
- ✓ Home Care

Finalé BRP3 Padded Slip-Style Bed Rail Protectors



ADVANTIFLEX PRESSURE CARE COVER

ADVANTIFLEX PU1100 POLYURETHANE PRESSURE CARE MATERIAL:

- ✓ Vapour permeable and breathable
- ✓ Anti-shearing
- ✓ Fire retardant, resistant to body fluids and bacteria
- ✓ Wipe-clean
- ✓ Heavy Duty PVC Inner Material – for maximum durability
- ✓ Optional: heavy duty PVC material for maximum durability

CONSTRUCTION:

- ✓ Closure style: velcro or press - for quick and simple fitment.
- ✓ Robust seams

25MM PADDED PROTECTIVE FOAM CORE


- ✓ High quality, open cell foam
- ✓ Anti-microbial treatment
- ✓ Padded on both sides of rail for maximum patient protection

APPLICATIONS

- ✓ Hospitals
- ✓ Aged Care
- ✓ Group Homes
- ✓ Home Care

Finalé Medical Bed Bolsters

Finalé Medical Bed Bolsters are customisable and available in a wide range of options to meet with varied patient needs. Available in core structures from a single layer foam to multi-layered, high care configurations, Finalé Bolsters are handmade from quality raw materials that meet and exceed clinical care requirements. Cover options include Heavy Duty PVC, Advantiflex Pressure Care Material and Mental Health Grade Material.

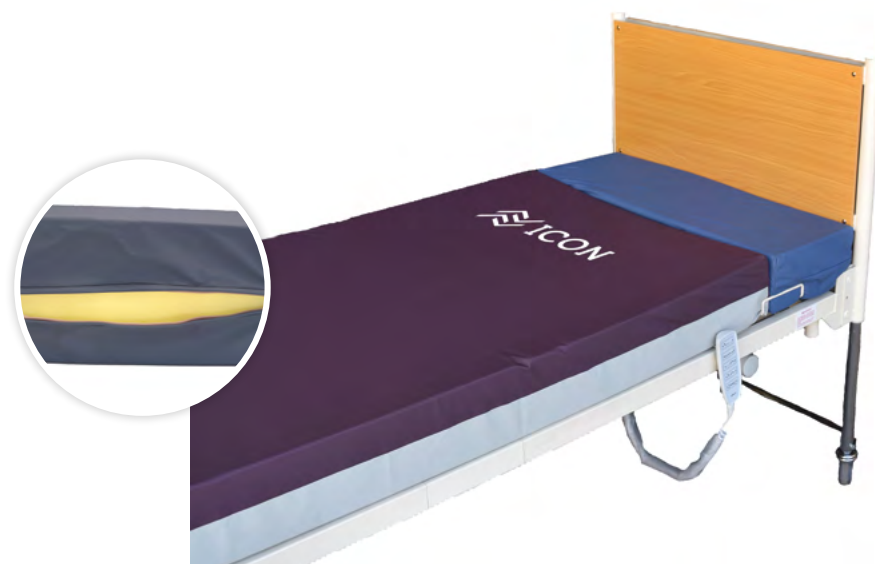

ADVANTIFLEX PRESSURE CARE COVER

ADVANTIFLEX PU1100 POLYURETHANE PRESSURE CARE MATERIAL:

- ✓ Vapour permeable and breathable
- ✓ Anti-shearing
- ✓ Fire retardant, resistant to body fluids and bacteria
- ✓ Wipe-clean
- ✓ Heavy Duty PVC Inner Material – for maximum durability
- ✓ Optional: heavy duty PVC material for maximum durability

CONSTRUCTION:

- ✓ Integrated robust zipper for cover removal
- ✓ Fluid penetration protection flap over zipper
- ✓ Robust seams

QUALITY FOAM CORE

- ✓ Supportive, durable
- ✓ Anti-microbial treatment

APPLICATIONS



- ✓ Mattress Extensions
- ✓ Positioning Aids

Finalé Medical Bed Bolsters

COVER	MATERIAL	SIZE (MM)	CODE: 150MM THICK	CODE: 180MM THICK	CODE: 200MM THICK
Sewn Seams	Advantiflex	880x200	BOL1S1015	BOL1S1018	BOL1S1020
		1050x200	BOL1KS1015	BOL1KS1018	BOL1KS1020
	Heavy Duty PVC	880x200	BOL3S1015	BOL3S1018	BOL3S1020
		1050x200	BOL3KS1015	BOL3KS1018	BOL3KS1020
HF Welded Seams	Advantiflex	880x200	BOL11S1015	BOL11S1018	BOL11S1020
		1050x200	BOL11KS1015	BOL11KS1018	BOL11KS1020
	Heavy Duty PVC	880x200	BOL13S1015	BOL13S1018	BOL13S1020
		1050x200	BOL13KS1015	BOL13KS1018	BOL13KS1020

Finalé Fall Protection Floor Mat Standard

Forté Fall Protection Floor Mats provide impact protection for patients at risk of falling from the bed surface. Care staff can have peace of mind that the Finalé Fall Mat utilises an extremely high impact absorbing foam core for superior patient protection, whilst meeting any infection control and WH&S concerns. Convenient, lightweight, handmade and durable, the Finalé Fall Mats are an essential part of any facility falls-management program.

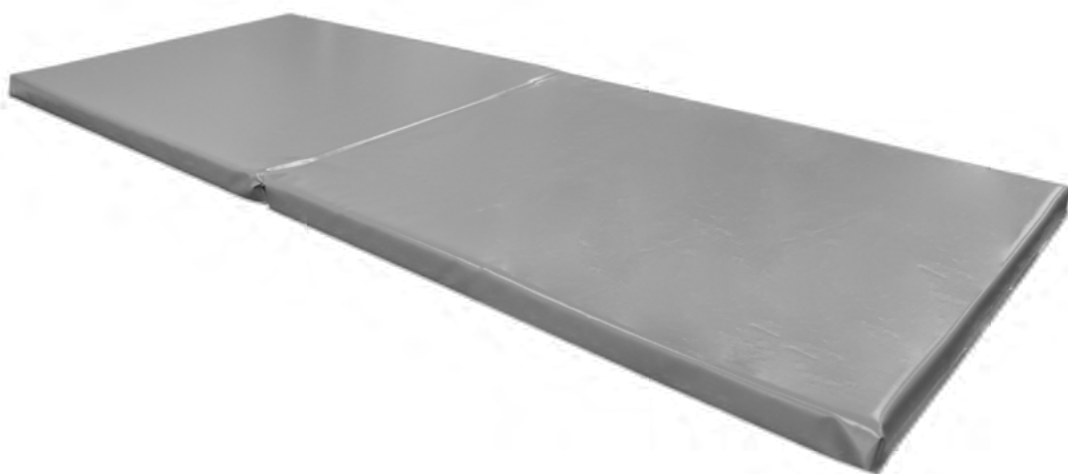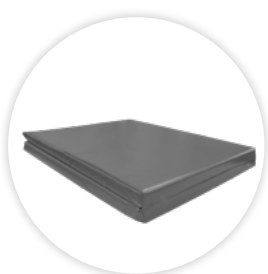

FULLY SEALED & ENCLOSED HEAVY DUTY PVC2 COVER

HEAVY DUTY, DURABLE PVC1120 MATERIAL:

- ✓ Wipe-clean, Hygienic
- ✓ 100% Waterproof
- ✓ Heavy Duty, Robust & Hardwearing
- ✓ 100% waterproof
- ✓ Anti-Microbial Protection
- ✓ Meets Source 5 Fire Retardant Standardsbody fluids
- ✓ Superior resistance to disinfection chemicals

CONSTRUCTION:

- ✓ Fully sealed and enclosed, no point of entrance to foam core
- ✓ Folds in half for storage
- ✓ Robust seams, Reinforced Critical Stress Points

DUAL DENSITY, HIGH IMPACT ABSORBING FOAM CORE

- ✓ Softens patients impact in event of a fall
- ✓ High quality, durable foam
- ✓ Folds for Storage

FOLDS FOR STORAGE

- ✓ Easy to relocate
- ✓ Compact and space-saving

Finalé Fall Protection Floor Mat - Standard

SIZE (MM)	CODE
1900x750x50	FMS1000

APPLICATIONS

- ✓ Hospitals
- ✓ Aged Care
- ✓ Home Care



Finalé Beveled Slimline Fall Protection Floor Mat

The Bevelled Slimline Finalé Fall Protection has been designed to provide ultimate protection for patients at risk of falling from the bed surface without jeopardising on WH&S and Infection Control. This advanced new technology is an essential part of any facilities falls management program.

FEATURES

Overall Size: 1800x700x40mm

Core: Dual Density, Impact Absorbing, Slim-Line Foam

- High Density Soft Interface
- Supportive EVA Foam Foundation
- Tri-fold Design - easy to carry & store
- Lightweight, great for manual handling concerns
- Beveled, to prevent trip hazard risks
- Slimline, low-profile with a re-inforced flat lay design - 40mm at maximum.
- Size meets the following relevant standards: VA 2014, Bowers et al 2008

Cover: Durable PVC3 Cover w/Non-Slip Patch

- Durable, Long Lasting PVC Coated Polyester
 - 630gsm
 - Provides Core Protection
 - Ultimate Durability
 - 100% waterproof to prevent ingress of fluid and infection
 - Fire retardant
 - Matte finish
 - Superior resistance to cleaning and disinfection chemicals
- High Performance Non-Slip/Anti-Skid Polyurethane Patch
 - Positioned on the underside of the mattress, Grips to floor surface
 - Vapour permeable/breathable
 - Resistant to body fluids, stain resistant
 - Wipe-clean, hygienic
 - Improves safety of care staff, Prevents mat from slipping/shifting
- Reflective Tape on all edges, to prevent trip hazard risks in poor lighting
- All Seams High Frequency Sealed and Welded
 - Electromagnetically fused seams means fluid cannot penetrate or soak through. The nature of the application of fall mats (used on the floor) is such that it is liable exposure to fluid spills and infection hazards
 - Ensuring spills and fluids can be completely wiped off and without lingering residue is important to ensuring a contamination free environment
 - Meets and exceeds infection control requirements
 - Reinforced with additional heavy duty stitching at stress points
 - No exposed sewn seams
- Fully Sealed & Enclosed
- Available in a wide variety of fabric colours - helps to distinguish between wards and ensure mats are not lost within facilities

Australian Made by Forté Healthcare

Images

