



ICON 135

PREMIUM 4-LAYER MEDIUM RISK
PRESSURE REDISTRIBUTION MATTRESS



The Icon Medium Risk Pressure Care Mattress Range provides true leadership in static pressure care in a preventative setting. Tried and trusted in facilities nation-wide, the Icon range combines outstanding durability, a high standard and consistency of care and excellent value. An extensive variety of configurations are available providing for a wide range of patient needs and requirements.

Icon Medical Mattress Characteristics




Airoform Immersion Interface

- » Envelops 'bony prominences'
- » Increases contact between surface and patient, dispersing and equalizing contact pressure
- » Provides excellent conformity to individual patient body profiles



Anatomical Support Region

- » Provides additional immersion and envelopment
- » Provides postural support for patient comfort
- » Resilient to optimize function, mobility and repositionings

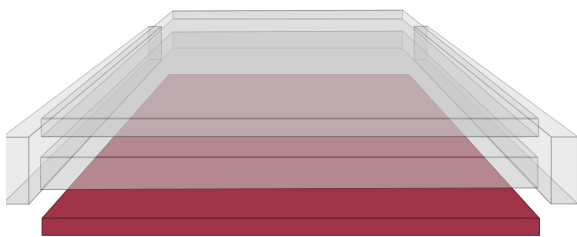

Firm Foundation and Base Support Region

- » Protects patient from 'bottoming-out' whilst in seated/semifowler/fowler position
- » Provides enhanced support for repositioning and transfer
- » Assists in primary care



Integrated Strengthened Sides

- » Provide lateral stability and reduce falls risk
- » Support patient transfer and ingress/egress
- » Assist primary care, provides caregivers a firm edge for support for repositioning
- » Protects patient from edge of bed surface



Excellent Immersion



Anatomical Support




Firm Foundation



Falls Protection

ICON I35 - 4-Layered Medium Risk Pressure Care Mattress

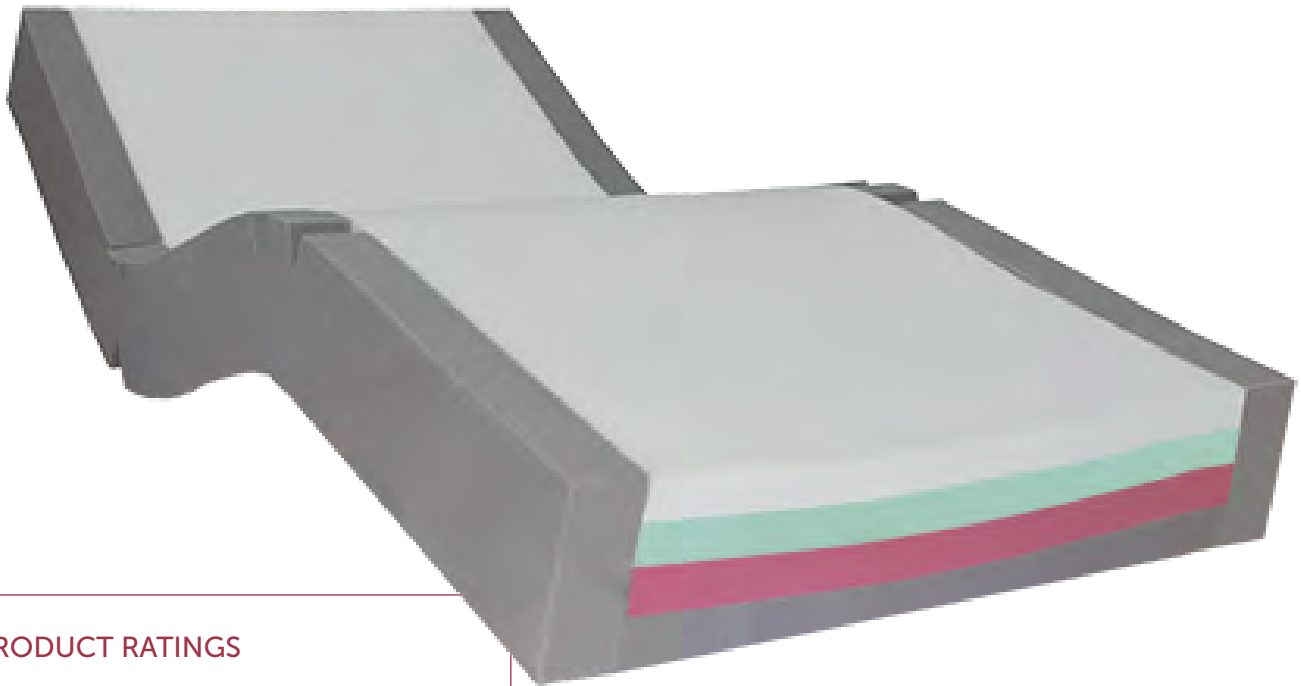

PRODUCT RATINGS

1 2 3 4 5	Shear & Friction
1 2 3 4 5	Repositioning
1 2 3 4 5	Infection Control
1 2 3 4 5	Firmness
1 2 3 4 5	Moisture & Microclimate
1 2 3 4 5	Performance & Longevity
1 2 3 4 5	Mobility & Function
MEDIUM RISK	Pressure Care Rating

APPLICATIONS

Hospital Wards
Aged Care
Palliative Care
Homecare

WEIGHT RANGE

40-200kgs

WARRANTY

5 Year Foam Core
2 Year Cover

CODE	DIMENSIONS (MM)	SIZE
I35S1000	1980x880x150	Single
I35SL1000	2030x880x150	Single, Long
I35KS1000	1980x1050x150	King Single
I35KSL1000	2030x1050x150	King Single, Long
I35EW1000	1980x1150x150	Extra-Wide
I35EWL1000	2030x1150x150	Extra-Wide, Long
I35D1000	1980x1350x150	Double
I35DL1000	2030x1350x150	Double, Long
I35Q1000	1980x1530x150	Queen
I35QL1000	2030x1530x150	Queen, Long
I35K1000	1980x1850x150	King
I35KL1000	2030x1850x150	King, Long
	Customisable	CAN BE MADE TO SUIT ANY BED




Airoform Memory breathe Immersion Interface Layer

- » High Temperature Stability to deliver consistent immersion, performance and comfort
- » Maximises Patient contact with Support Surfaces Providing Pressure Redistribution
- » Alleviates localised pressure points
- » Depth calibrated to maximise functionality and patient mobility
- » Unique Airoform Breathe ventilation system minimises heat entrapment and enhances core circulation providing a cooler sleep surface
- » Visco-elastic/Memory Foam cell structure provides low-shear surface
- » Firmness: Soft



4-Layered, Graded Pressure Relieving Core

- » High specification foam for longevity and performance
- » Provides graded postural support and patient comfort
- » Anti-microbial treatment applied to foam maximises longevity and core micro-climate
- » Excellent core stability to encourage centralised patient positioning to minimise falls risk, entrapment and enhance pressure relief

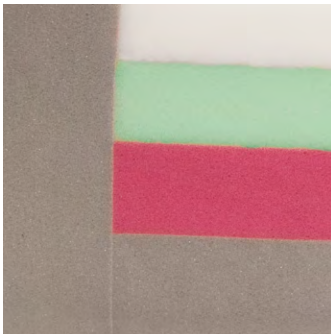


Premium Firm foundation Layer

- » Provides protection against 'bottoming-out' during transfer, repositioning and whilst in the inclined/fowler position
- » Provides exceptional comfort, support for repositioning and patient care



hinged strengthened sides


- » Increase lateral stability of support surface
- » Improves patient safety and encourages central patient positioning
- » Make patient transfers easier
- » Unique hinge profile reduces foam fatigue and allow mattress to easily conform to varied bed positions
- » Enables ease of primary patient care

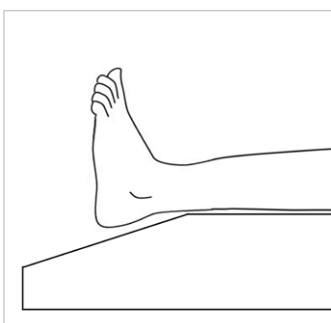

Sloped heel pressure transfer zone

- » Enhances pressure relief & protection for critical at-risk heel region
- » Offloads pressure from the patients heels and transfers it to calf area
- » Unobtrusive, does not require any adjustment

Pressure Mapping

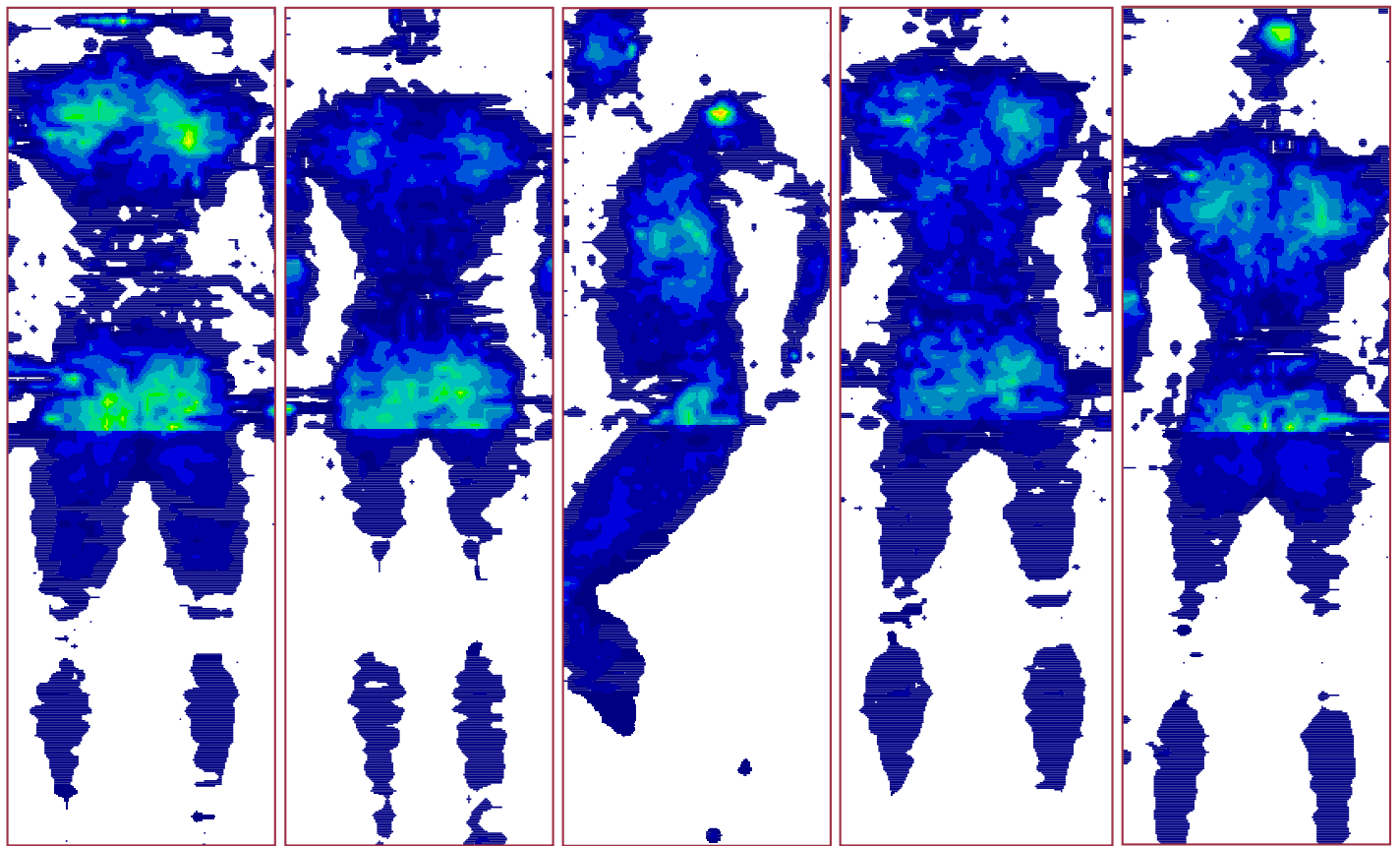


FOWLER
- ELEVATED LEGS

FOWLER

LATERAL RECUMBANT

SUPINE - ELEVATED LEGS

SUPINE

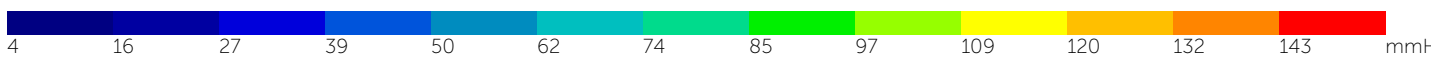

NOTES - Pressure Mapping Parameters

Bed Surface: Alrick, 5 Function Electric Hospital/Homecare Bed

Ambient Temperature: 21°C

Duration: Pressure Mapping was captured after 6 hours of Subject Laying Upon Support Surface

Subject: 98kg, Male, Height: 188cm

Tekscan Pressure Measurement System 7.20C
SENSOR_TYPE 5400D
ROWS 104
COLS 34
SENSEL_AREA 2.89008 cm²
SECONDS_PER_FRAME 0.010319
MICRO_SECOND 0
UNITS mmHg

Covers



New! Luxury Quilted Cotton Cover

- » Available on all Forté Mattresses
- » Contact Forté for further information

Mattress & Cushion Covers

	ADVANTIFLEX CLASSIC	PREMIFLEX CLASSIC	PREMIFLEX ULTRA
High Frequency Welded Seams	YES	YES	YES
Integrated Zipper	YES	YES	YES
Waterfall Zipper Protection Flap	YES	YES	YES
4-Way Stretch Patient Contact Surface	YES	YES	YES
Cover Warranty	2 Years	4 Years	4 Years
Chlorine / Sodium Hypochlorite Resistance	1,000PPM	10,000PPM	10,000PPM
Heavy Duty PVC1100 Base	YES		
Heavy Duty PVC1200 Base		YES	
Heavy Duty DoubleCoated Polyurethane Base			YES
Premium Comfort, Non-Abrasive and Durable Upper	YES		
Superior Durability, Increased Shear Protection, Ultimate Comfort Layer		YES	YES
Full Vapour Permeability	YES	YES	YES

Cover Features

Advantiflex Classic Pressure Care Cover

Belgium made polyurethane coated upper and heavy duty PVC Base ensures that Forté Advantiflex Covers deliver consistent care, comfort and durability. Utilising premium quality medical grade fabrics, Forté Combination Covers are a tried and proven medical grade Pressure Relieving Cover.

Premiflex Classic Pressure Care Cover

Forté Premiflex Classic Covers provide outstanding infection control, hygiene and high performance within the Acute Care Setting. Manufactured with a High-Quality Belgium Halogen Free Polyurethane 4-Way Stretch, Moisture Vapour Permeable Anti-Shearing Interface Material, the Premiflex Cover is an effective Pressure Injury risk mitigation mattress cover.

Premiflex Ultra Pressure Care Cover

Forté Premiflex Ultra Covers provide the ultimate protection and performance for mattresses in an acute care setting. Superior pressure injury protection is achieved through 4-way stretch non-hammocking and anti-shearing characteristics, Premiflex Ultra Covers have a durable double coated polyurethane base that is moisture vapour permeable allowing foam longevity by allowing the core to breathe.




High Frequency Sealed and Welded Seams

- » Electromagnetically fused seams protecting against risk of penetration of fluid
- » Meets and exceeds infection control requirements
- » Reinforced with additional heavy duty stitching at stress points
- » 100% sealed to prevent ingress of fluid or infection



Integrated Zipper

- » Allows easy cover removal
- » Robust and durable
- » High frequency welded attachment for optimal protection against fluid ingress

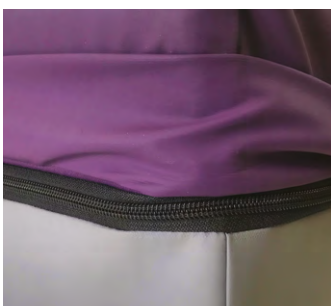

Waterfall Zipper Protection Flap

- » Protects against fluid ingress at point of zipper
- » Meets infection control guidelines



GLOBAL TECHNOLOGY
AUSTRALIAN MADE

All mattresses are made by
Forté Healthcare in Armidale
NSW in accordance with all
Australian manufacturing
and clinical standards.